


Sure Status COVID-19 IgG/IgM Card Test

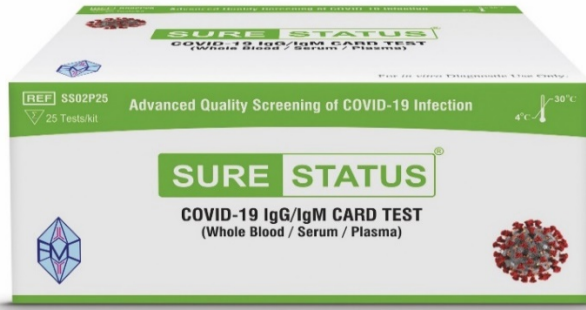

Sure Status COVID-19 IgG/IgM Card Test is a rapid immunochromatography test designed for the qualitative detection of specific IgM and IgG to SARS-CoV-2 in human Blood/Serum/Plasma.

- Rapid testing for SARS-CoV-2 antibodies within 15-20 minutes
- Add 2 drops (20 μ L) of whole blood or 1 drop (10 μ L) of serum/ plasma to the specimen well
- Suitable for Point of Care Testing.
- No need for extra equipment, Easy to use
- Specimen : Whole blood/Serum/Plasma
- Based on Colloidal Gold immunochromatography (Visually read lateral flow immunoassay, Stable at 4-30°C)
- Shelf Life: 24 months
- Package: 25 Tests / kit
- Highest Sensitivity and Specificity

Ordering Information

Cat. No.	Product	Storage temperature	Pack size
SS02P25	Sure Status COVID-19 IgG/IgM Card Test	4-30°C	25T/kit

Order at info@premiermedcorp.com